

Quick Start Guide




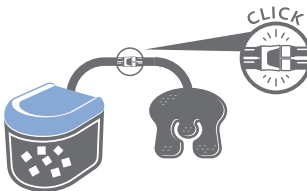
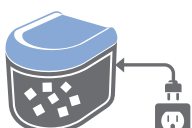
DONJOY®

DonJoy® IceMan CLASSIC^{3™} Cold Therapy Unit with Universal Wrap-On Pad - Knee Application




- THIS DEVICE CAN BE COLD ENOUGH TO CAUSE SERIOUS INJURY.
- DO NOT USE THIS DEVICE WITHOUT A PRESCRIPTION FROM A PHYSICIAN.
- READ AND UNDERSTAND ALL WARNINGS AND INSTRUCTIONS FOR USE BEFORE USING THIS DEVICE. ADDITIONAL WARNINGS APPEAR IN THE INSTRUCTIONS FOR USE.



Directions for use of the IceMan CLASSIC^{3™}

-  Add ice to fill line.
-  Add cold water to fill line.
-  Place lid on device making sure the label is facing up. Secure the lid by raising the handle, which will engage the lid locking mechanism.
-  Connect cold therapy hose to cooling pad. To ensure a reliable connection, "snap" or "click" hoses together into place so that the fit is tight and snug.
-  To turn device on, insert cord into connection on the back of the device and plug power supply into the wall outlet.

Directions for use of the Wrap-On Pad

-  Place barrier on knee then place wrap-on pad on knee with tube up as shown. Place round flap on top of patella.
-  Take strap and wrap around knee from one side of the pad over the patella. Wrap should be snug but not tight.
-  Secure long strap to itself to keep everything flat.

- Apply cold pad and barrier to affected area, making sure to use your prescription.
- **WARNING!** - When applying the cold pad, DO NOT let any part of the cold pad touch your skin.
- Always use with a barrier between your skin and the cold pad.
- Check for moisture on the barrier between your skin and cold pad.
- If moisture is present on the barrier, immediately discontinue use of this device.
- Do not use this device without a prescription from a physician.

Product Support | +1.888.405.3251

WARRANTY: DJO, LLC will repair or replace all or part of the unit and its accessories for material or workmanship defects for a period of six months from the date of sale.

Individual results may vary. Neither DJO Global, Inc. nor any of its subsidiaries dispense medical advice. The contents of this sheet do not constitute medical, legal, or any other type of professional advice. Rather, please consult your healthcare professional for information on the courses of treatment, if any, which may be appropriate for you.



T 800.336.6569 D 760.727.1280 F 800.936.6569
2900 Lake Vista Drive | Lewisville, TX 75067 | U.S.A.
DJOGlobal.com/donjoy

Copyright © 2020 by DJO, LLC
MKT00-9677 Rev A